


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Vildagliptin (racemic)

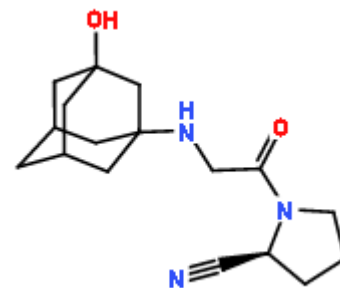

CHROMATOGRAM

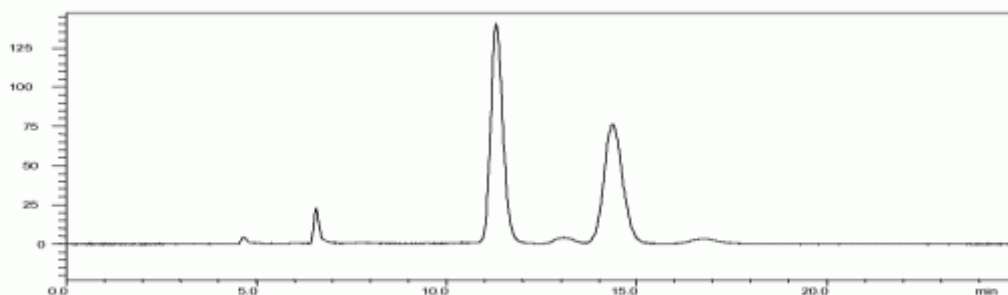

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	ethanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.3	0.89	1 - 2	4.02
2	14.3	1.39		