


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Kynurenine

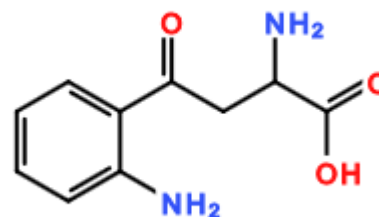

CHROMATOGRAM

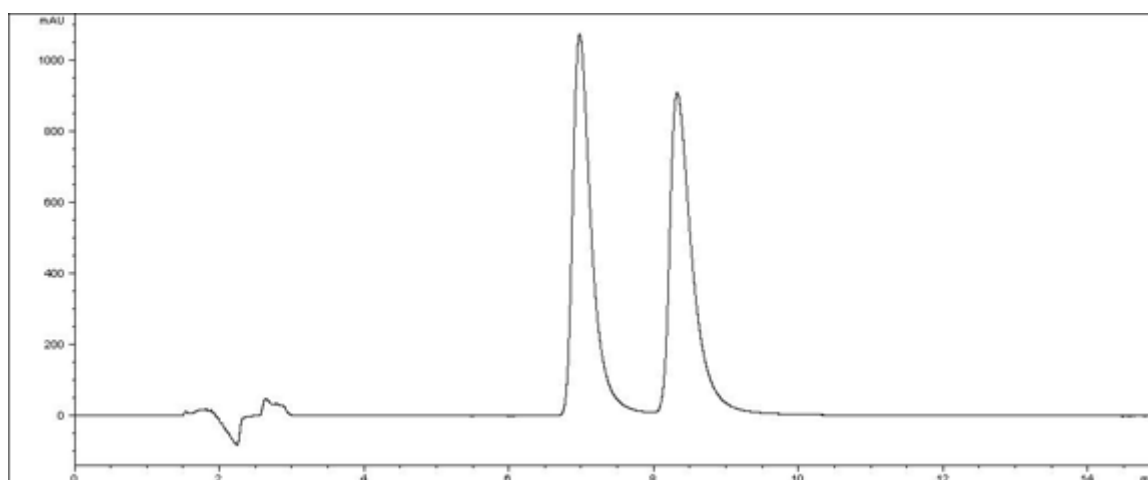

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	25mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order	L/D		

No.	Rt	K'	Alpha Value	Resolution
1	7.0	3.38	1 - 2 1.25	2.79
2	8.3	4.23		