


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Atenolol

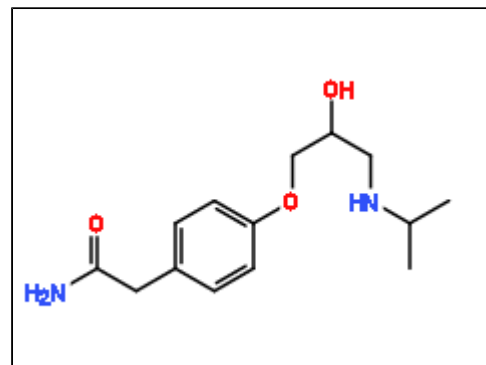

CHROMATOGRAM

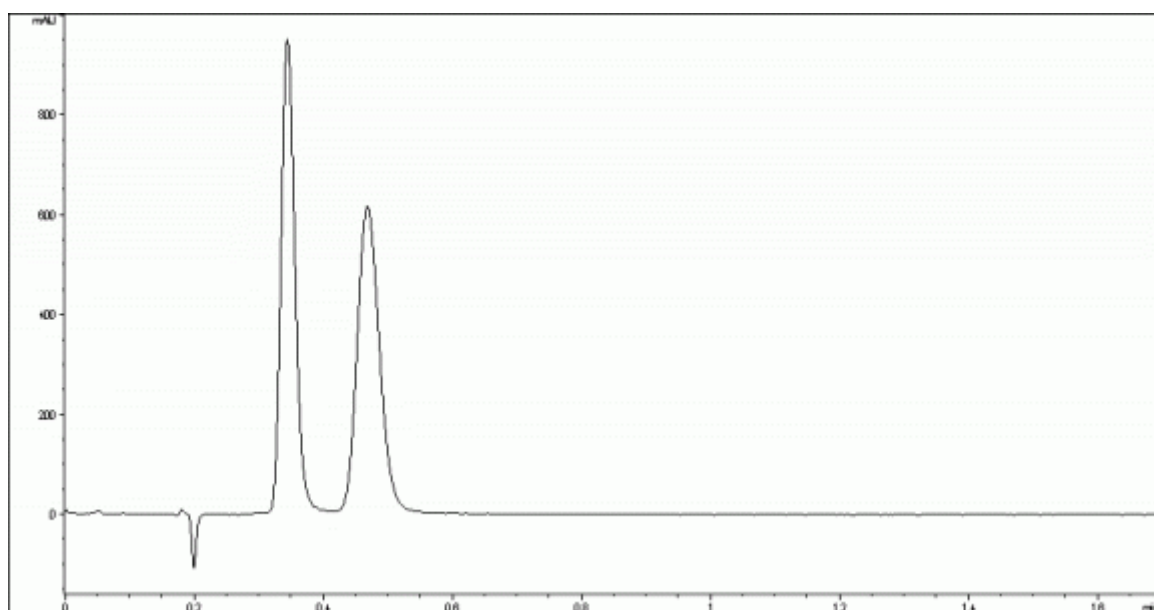

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	acetonitrile / 2-propanol / diethylamine = 95 / 5 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	4	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	1.27	1 - 2	2.47
2	0.5	2.09		