


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Dropropizine

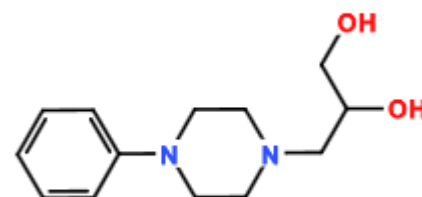

CHROMATOGRAM

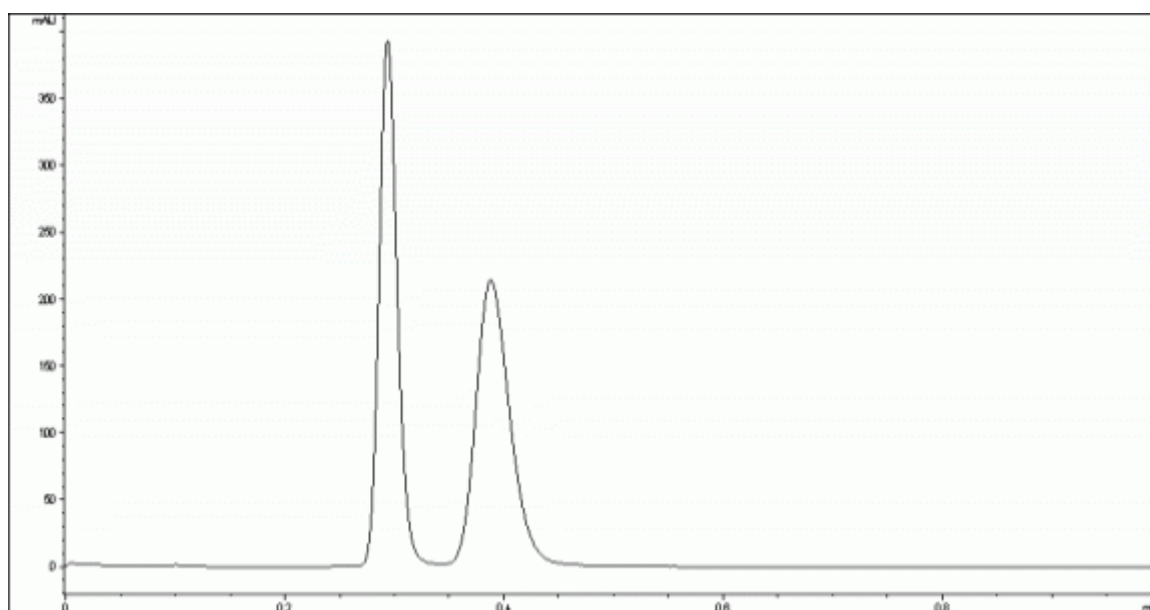

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	ethanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	3	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	0.46	1 - 2	2.02
2	0.4	0.92		2.21