


APPLICATION RESULTS




COMPOUND IDENTIFICATION

Compound Name
Ofloxacin

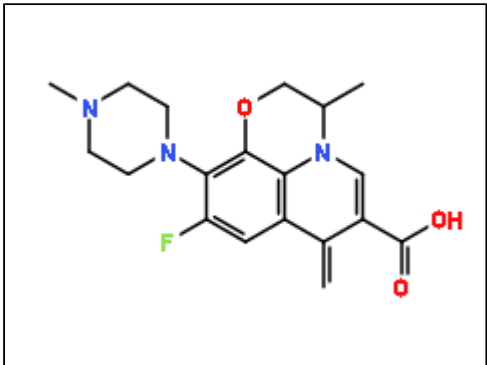

CHROMATOGRAM

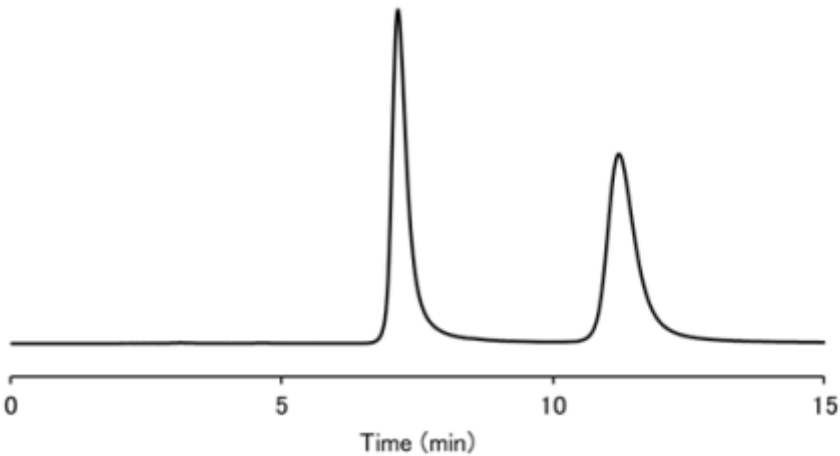

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OX-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	20mM ammonium acetate aq. / acetonitrile = 65 / 35		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	1	Sample Concentration	0.5 mg/mL
Detection	UV-VIS 291 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.1		1 - 2	5.68
2	11.2			