


APPLICATION RESULTS




COMPOUND IDENTIFICATION

Compound Name

5,5',6,6',7,7',8,8'-Octahydro-1,1'-bi-2-naphthol

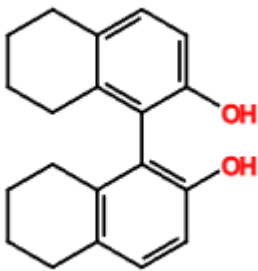

CHROMATOGRAM

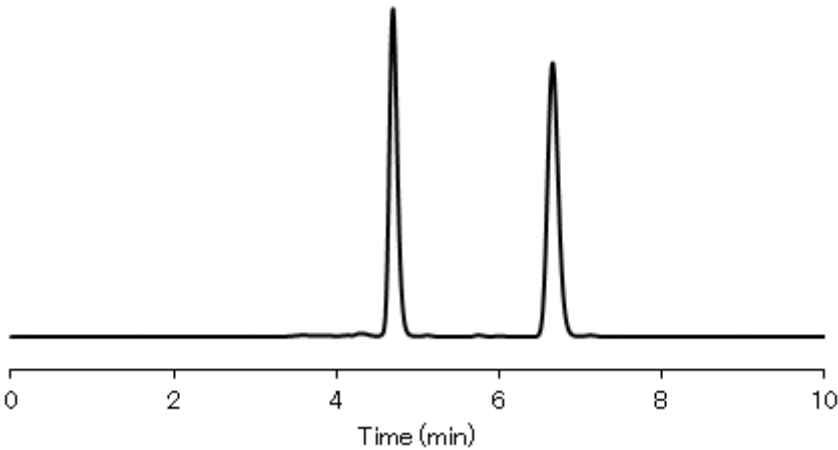

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 30 / 70		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	0.5 mg/mL
Detection	UV-VIS 285 nm		
Elution Order	第一ピーク : R体 第二ピーク : S体		

No.	Rt	K'	Alpha Value	Resolution
1	4.7		1 - 2	9.14
2	6.7			