


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

gamma-(4-Fluorophenyl)-gamma-butyrolactone

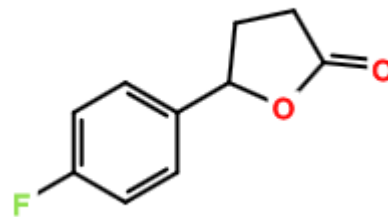

CHROMATOGRAM

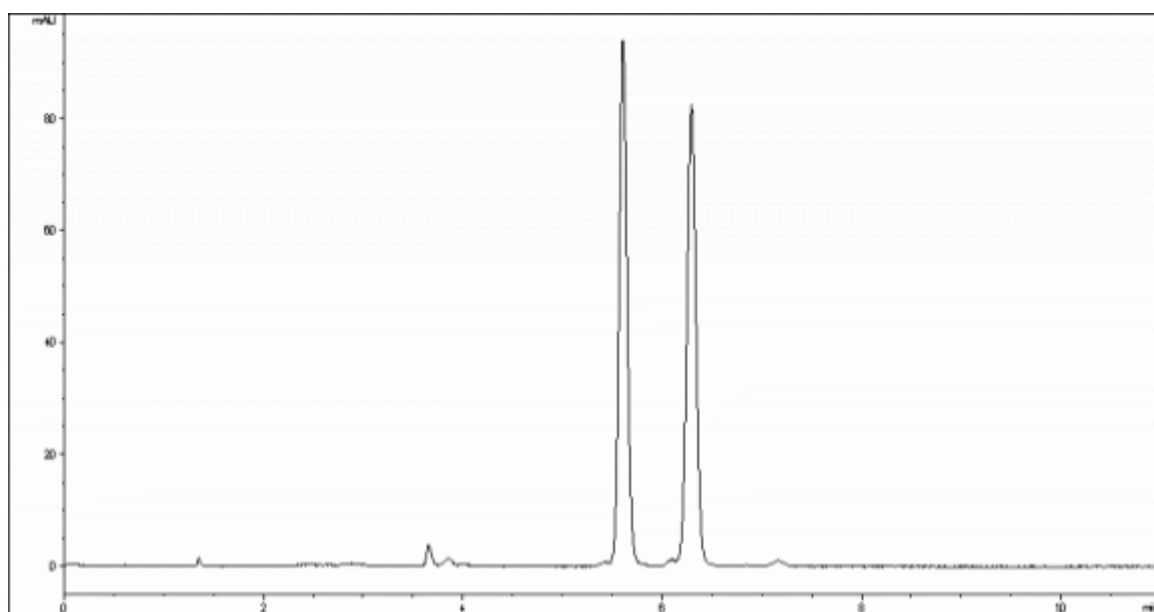

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] AY-3		
Column Size (Particle Size)	4.6 x 150 mm (3 μm)		
Mobile Phase Composition	ethanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.6	2.17	1 - 2 1.18	4.43
2	6.3	2.56		