


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Metalaxyl

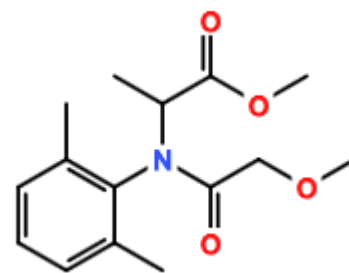

CHROMATOGRAM

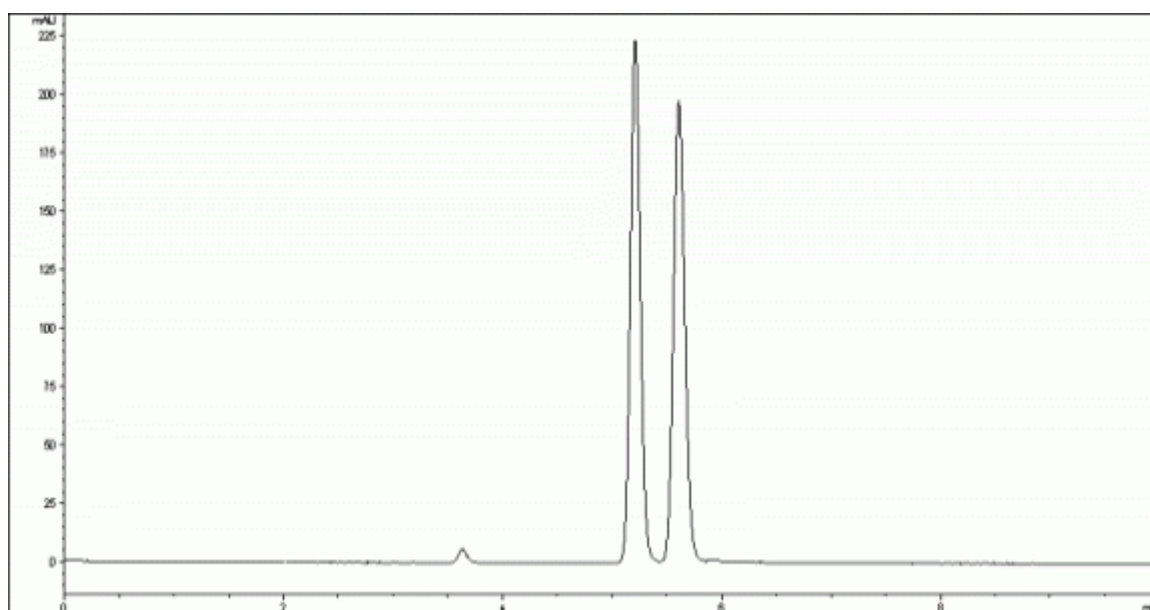

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OZ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	ethanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.2	2.08	1 - 2	2.26
2	5.6	2.32		