


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Labetalol

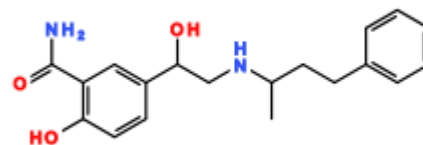

CHROMATOGRAM

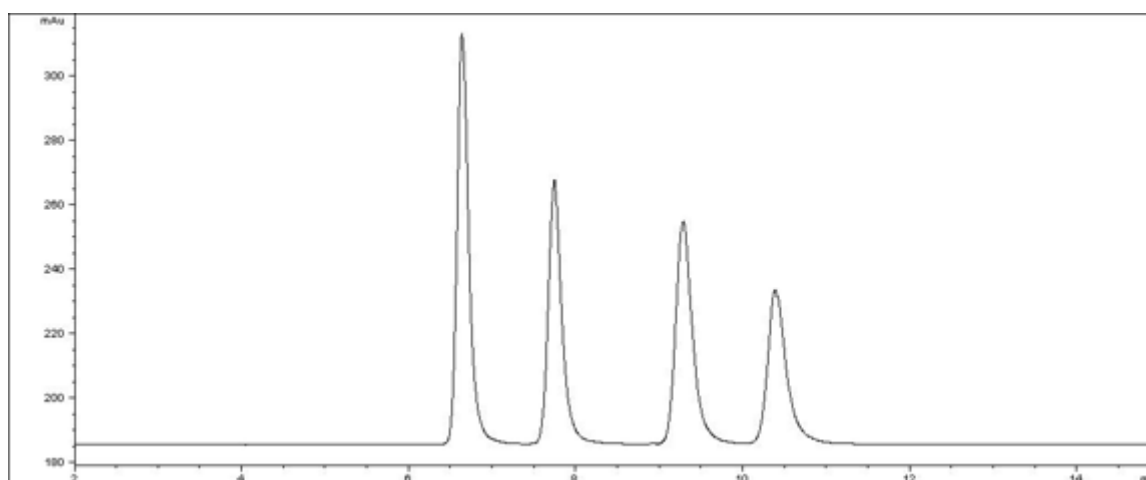

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	12.5mM triethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.6	2.88	1 - 2	4.22
2	7.7	3.52	2 - 3	4.88
3	9.3	4.42	3 - 4	2.99
4	10.4	5.06		