


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sertraline tetralone

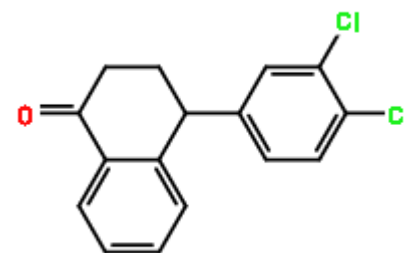

CHROMATOGRAM

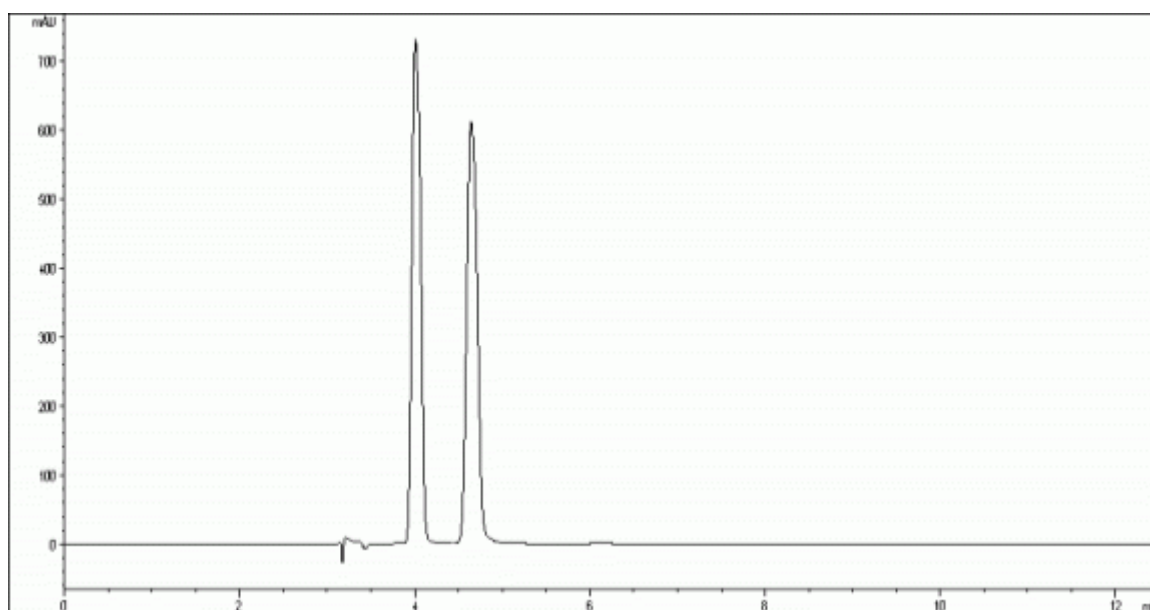

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AZ-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	acetonitrile / 2-propanol = 95 / 5		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.0	0.35	1 - 2	3.02
2	4.7	0.57		