


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Indapamide

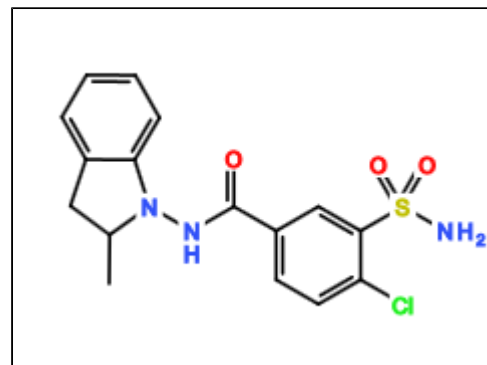

CHROMATOGRAM

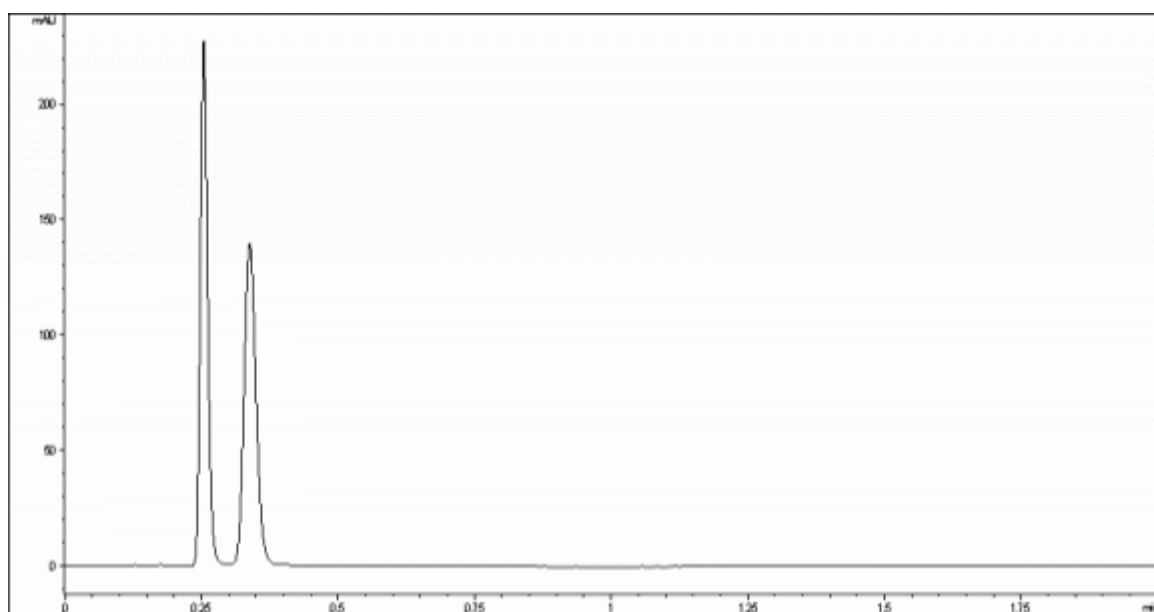

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OD-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	methanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	3	Temperature (°C)	40
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	1.14	1 - 2	2.65
2	0.3	1.85		