


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Isoamarine

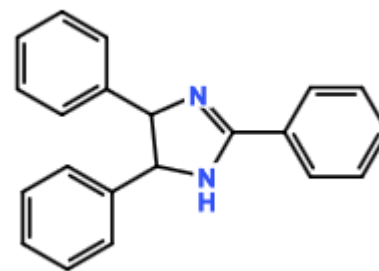

CHROMATOGRAM

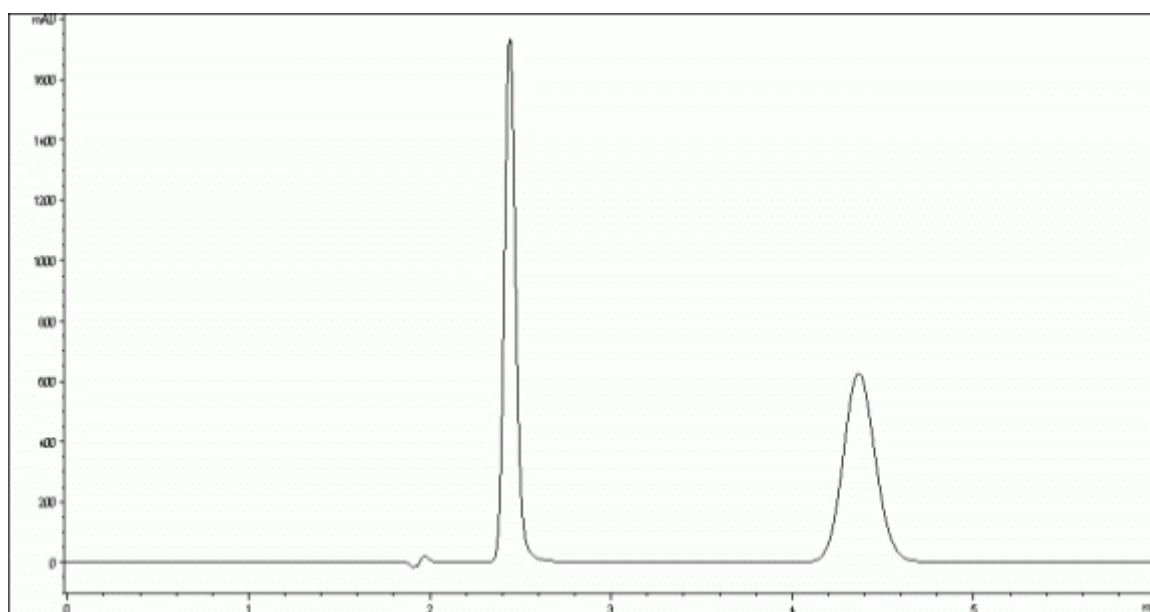

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OJ-H		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	methanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.4	0.31	1 - 2	8.84
2	4.4	1.34		