


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Valine

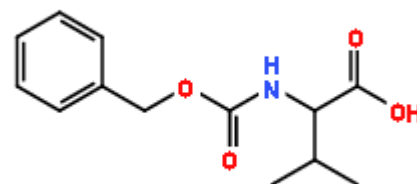

## CHROMATOGRAM

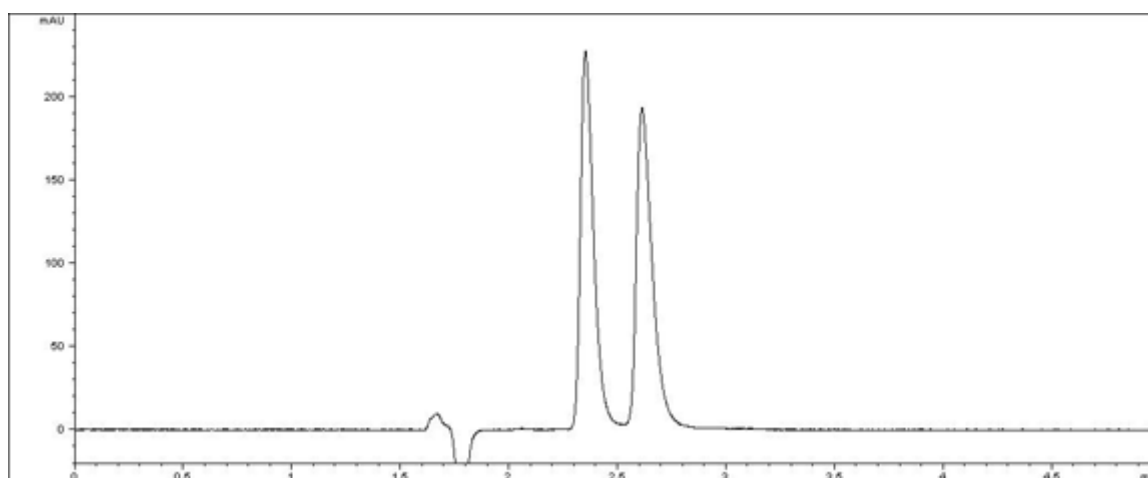

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	10
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 218 nm		
Elution Order	D/L		

No.	Rt	K'	Alpha Value	Resolution
1	2.4	0.41	1 - 2	2.19
2	2.6	0.56		