


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Methocarbamol

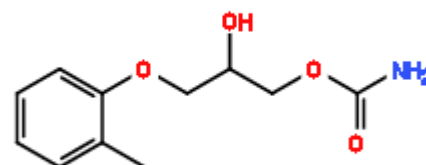

CHROMATOGRAM

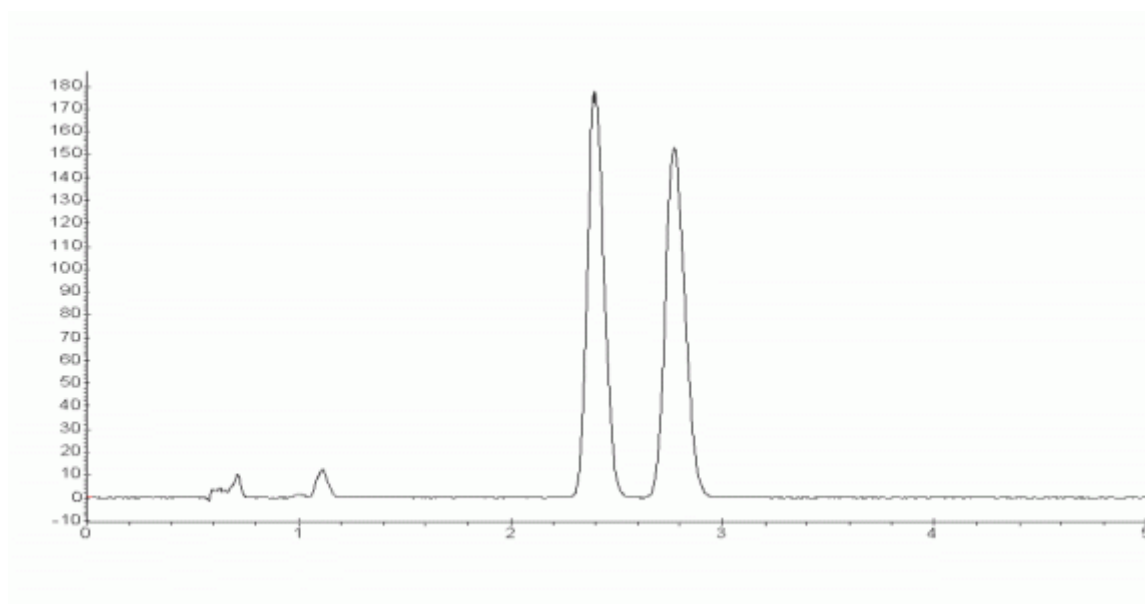

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / methanol = 80 / 20		
Chromatographic Mode	SFC		
Flow rate (mL/min)	5	Temperature (°C)	35
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.4	3.00	1 - 2	1.21
2	2.8	3.62		