


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-FMOC-DL-Phenylalanine

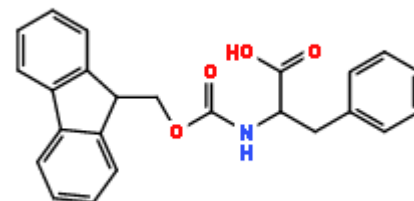

## CHROMATOGRAM

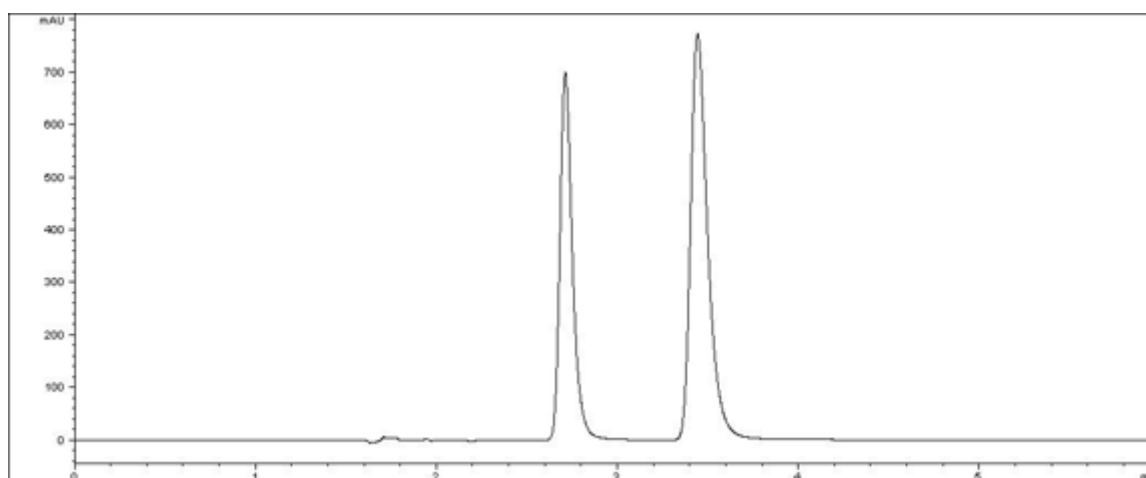

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order	D/L		

No.	Rt	K'	Alpha Value	Resolution
1	2.7	0.65	1 - 2	5.16
2	3.4	1.10		