


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Benzyl mandelate

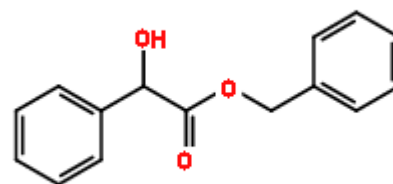

CHROMATOGRAM

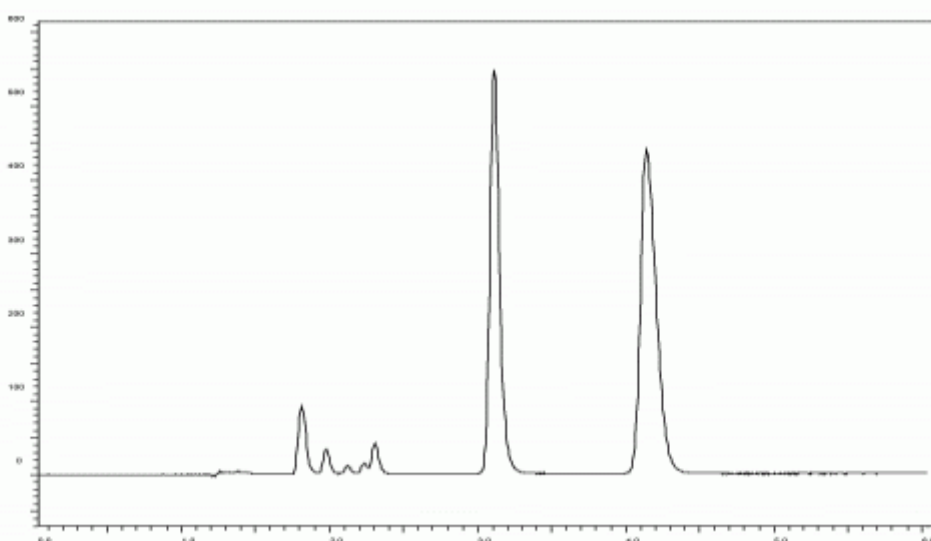

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] AY-H		
Column Size (Particle Size)	4.6 x 250 mm (5 μ m)		
Mobile Phase Composition	CO ₂ / methanol = 85 / 15		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3	Temperature ($^{\circ}$ C)	RT
Injection Volume (μ L)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.1	2.11	1 - 2	6.60
2	4.1	3.14		