


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Flavanone

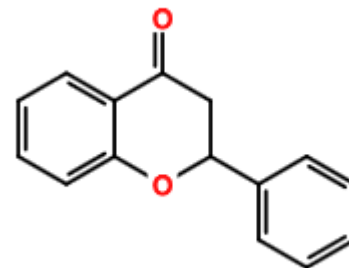

CHROMATOGRAM

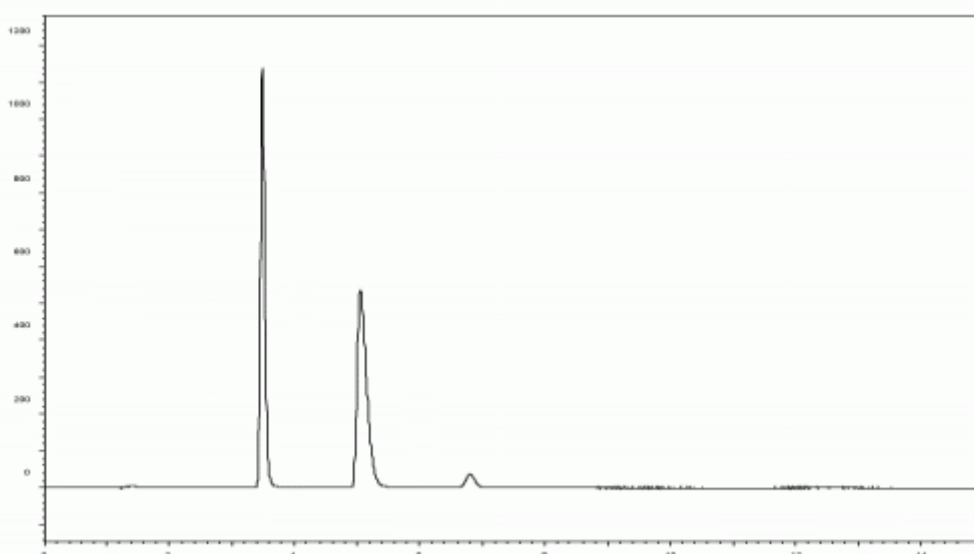

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] AY-H		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	CO ₂ / methanol = 85 / 15		
Chromatographic Mode	SFC		

Flow rate (mL/min)	3	Temperature (°C)	RT
Injection Volume (μL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	2.49	1 - 2	7.82
2	5.1	4.05		