


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Na,Ne -di-CBZ-DL-lysine

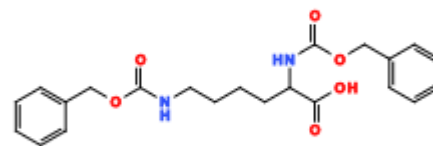

## CHROMATOGRAM

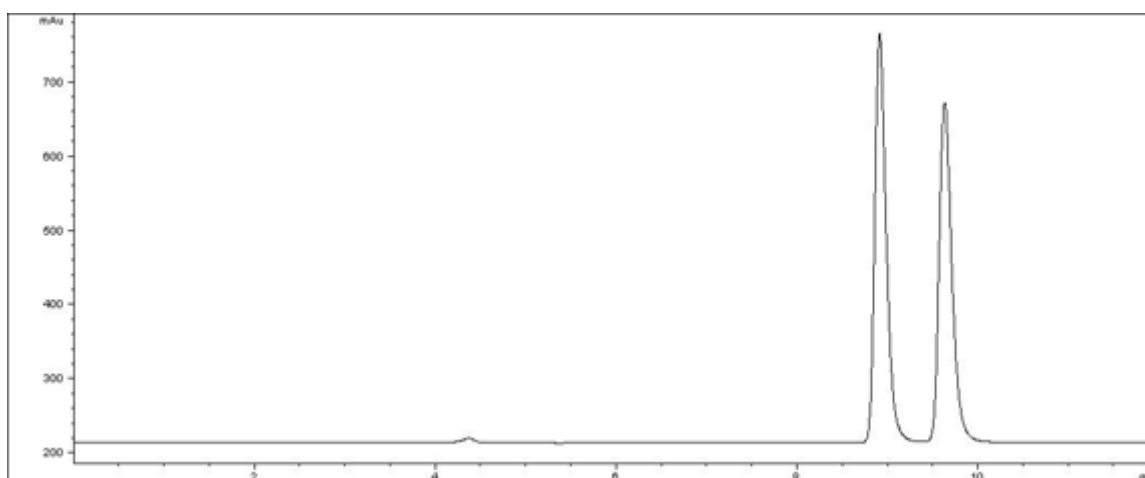

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 µm)		
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	10
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.9	2.14	1 - 2	3.06
2	9.6	2.39		