


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5-(4-Hydroxyphenyl)hydantoin

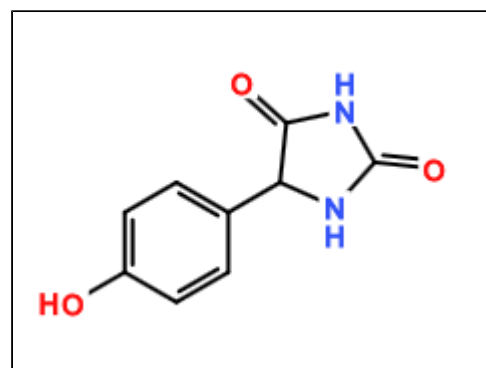

CHROMATOGRAM

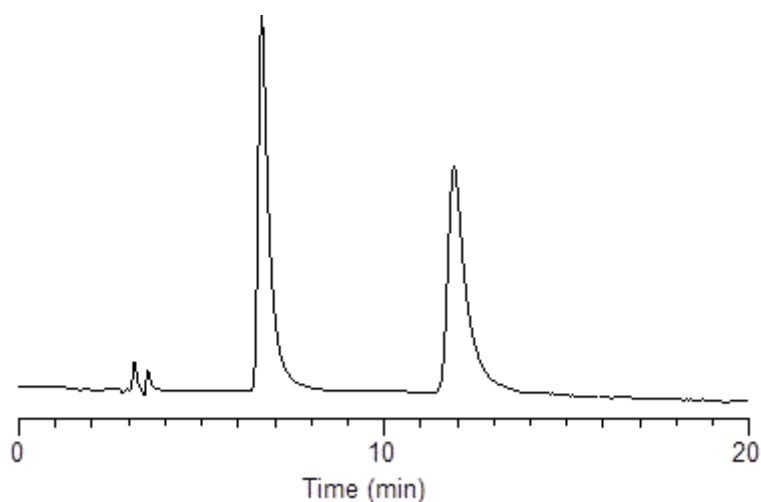

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol = 70 / 30		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.6	1.20	1 - 2	8.14
2	11.9	3.00		