


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sulindac methyl ester

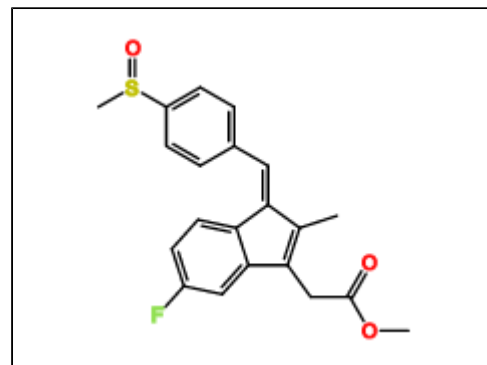

CHROMATOGRAM

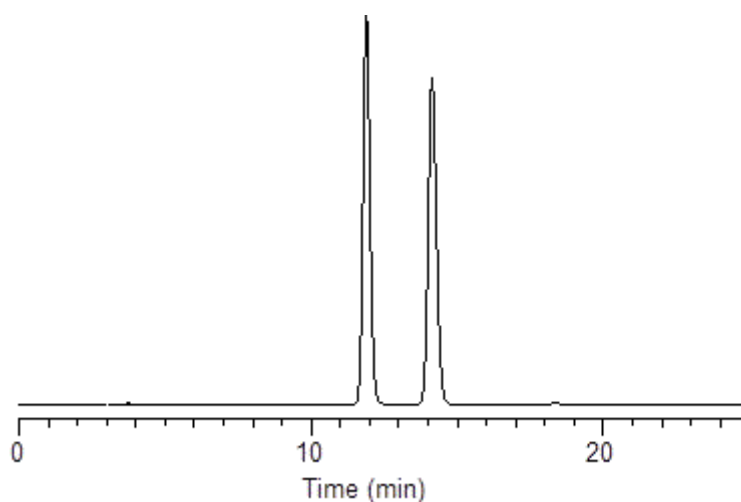

ANALYTICAL REPORT

| | | | |
|-----------------------------|------------------------------|----------------------|----|
| Column Trade Name | CHIRALPAK® IA | | |
| Column Size (Particle Size) | 4.6 x 250 mm (5 µm) | | |
| Mobile Phase Composition | n-hexane / ethanol = 80 / 20 | | |
| Chromatographic Mode | Normal Phase | | |
| Flow rate (mL/min) | 1.0 | Temperature (°C) | 25 |
| Injection Volume (µL) | 0.01 (mg) | Sample Concentration | |
| Detection | UV-VIS 254 nm | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|------|------|-------------|------------|
| 1 | 11.9 | 3.00 | 1 - 2 | 4.74 |
| 2 | 14.1 | 3.70 | | |