


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cetirizine

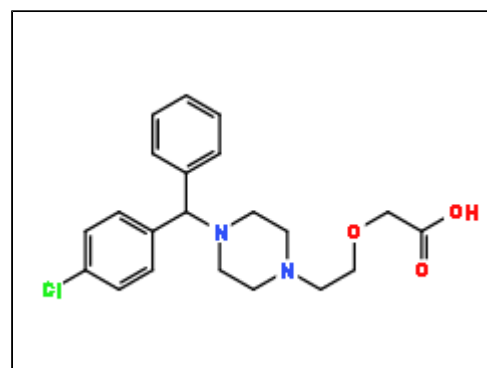

## CHROMATOGRAM

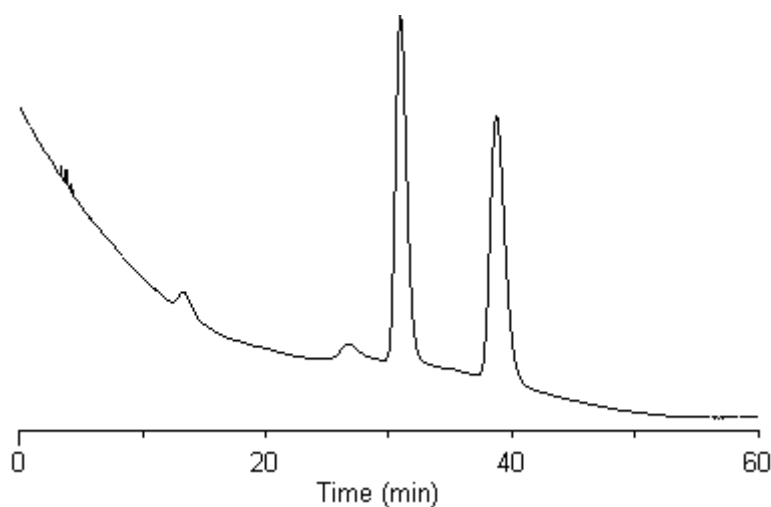

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IE 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / acetic acid / diethylamine = 80 / 20 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	30.9	9.30	1 - 2	4.00
2	38.7	11.90		