


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tiaprofenic acid

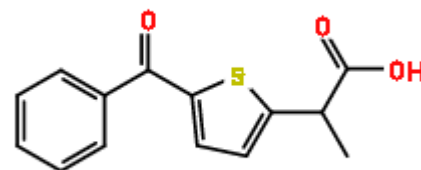

CHROMATOGRAM

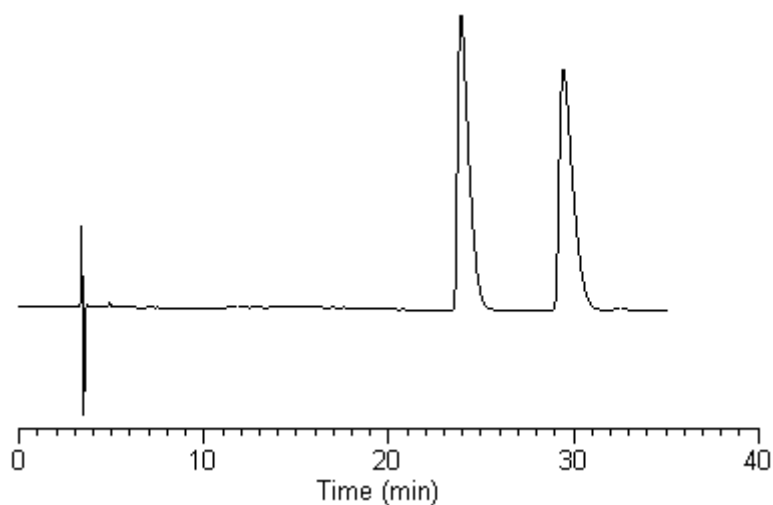

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IE 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	23.9	7.00	1 - 2 1.26	4.59
2	29.5	8.80		