


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Tryptophan

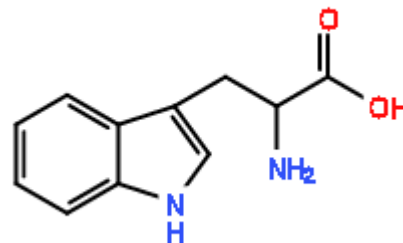

## CHROMATOGRAM

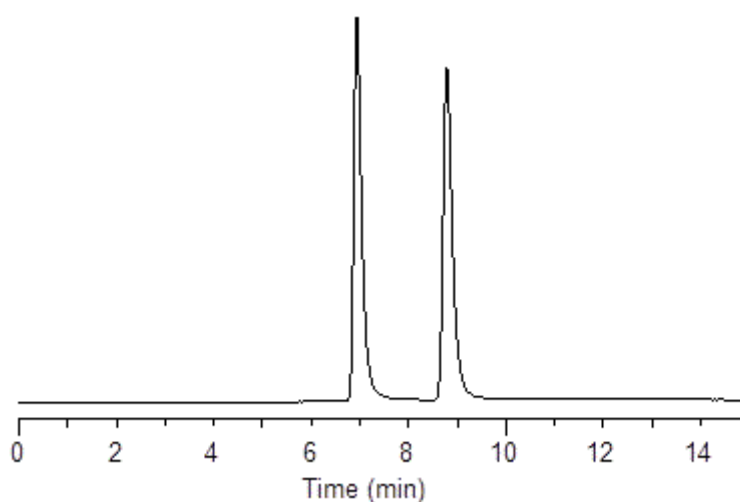

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 4 x 250 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.9	1.80	1 - 2	6.01
2	8.8	2.50		