


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

2-(4-Hydroxyphenoxy)propionic acid

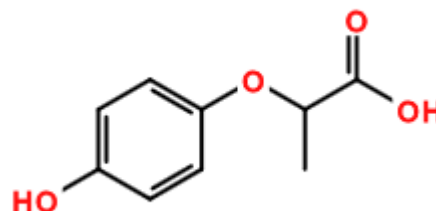

CHROMATOGRAM

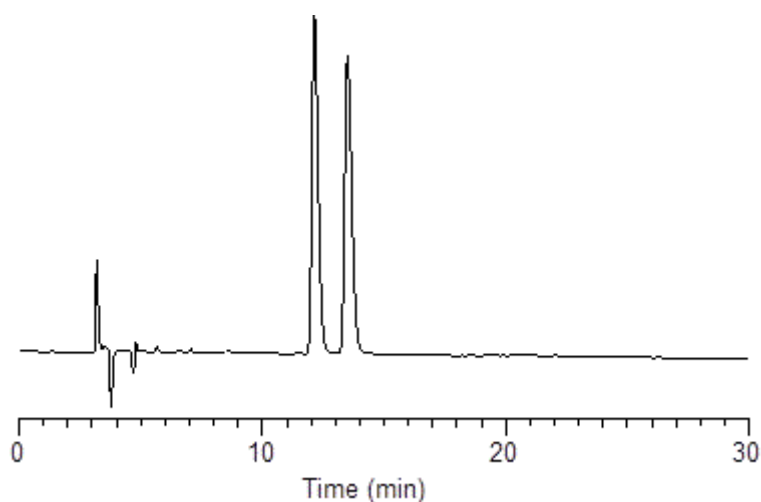

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 93 / 7 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.1	3.00	1 - 2	2.71
2	13.5	3.50		