


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

1,5-Dihydroxy-1,2,3,4-tetrahydronaphtalene

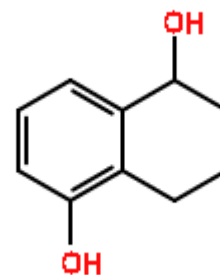

## CHROMATOGRAM

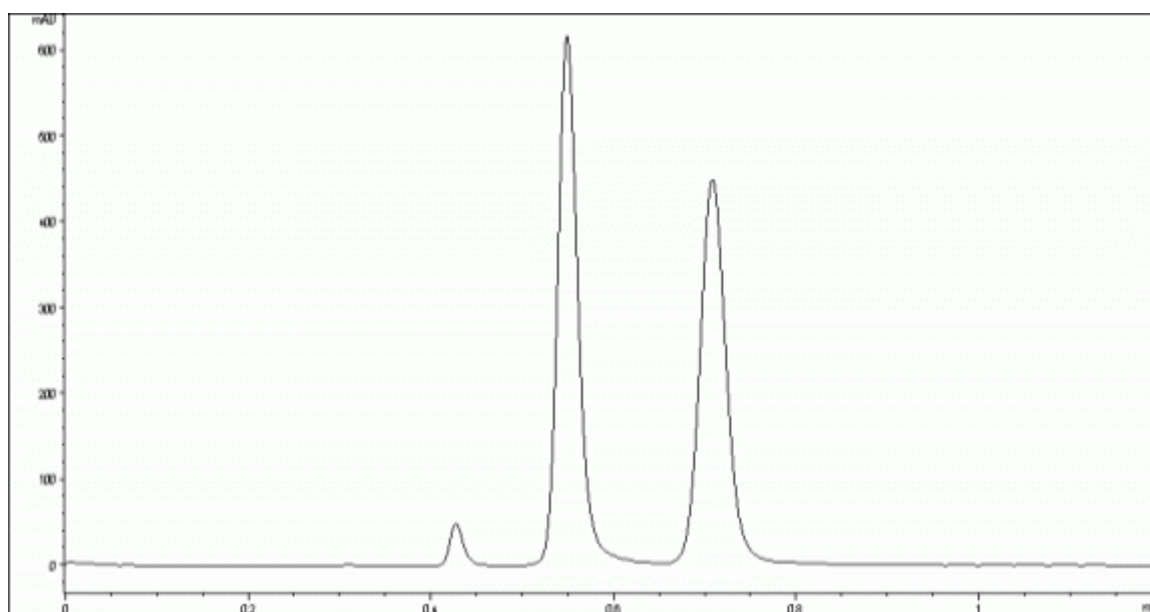

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	2	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.6	0.91	1 - 2	3.54
2	0.7	1.46		