


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

alpha-Acetyl-alpha-methyl-gamma-butyrolactone

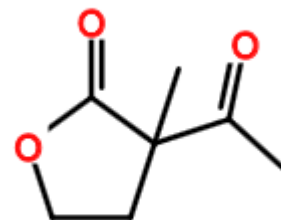

## CHROMATOGRAM

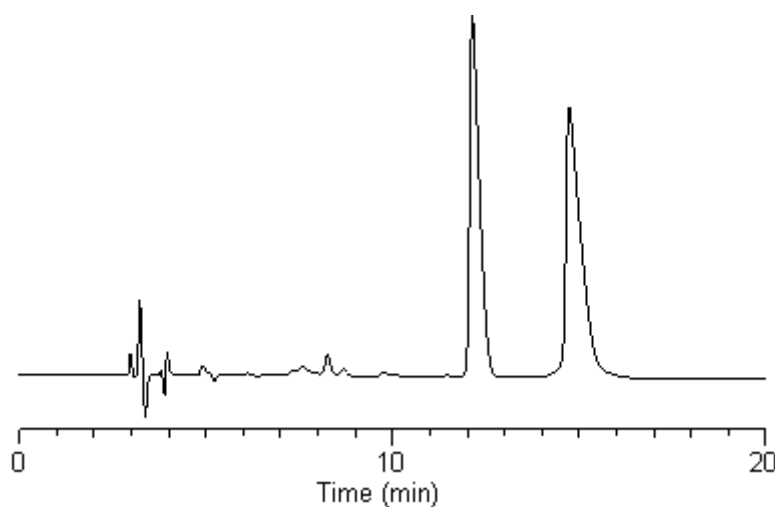

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol = 80 / 20		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.1	3.00	1 - 2	1.28
2	14.7	3.90		