


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Alvimopan

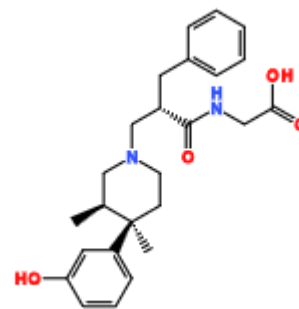

## CHROMATOGRAM

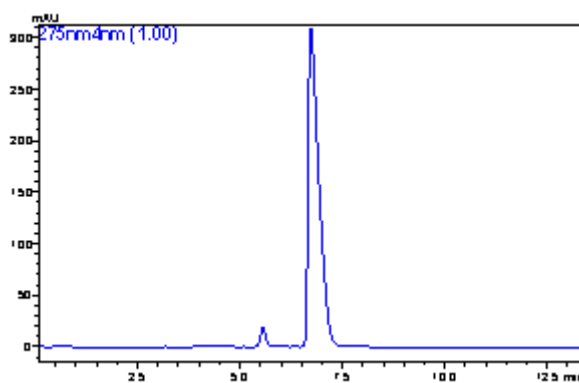

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / triethylamine = 75 / 25 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order	R, S		

No.	Rt	K'	Alpha Value	Resolution
1	5.6	0.85	1 - 2	1.46
2	6.7	1.24		3.12