


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Bepotastine Besilate

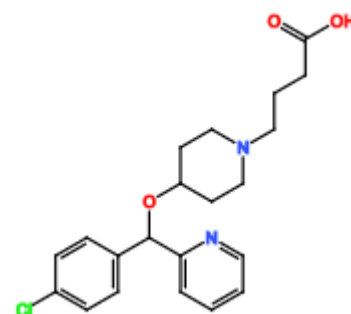

CHROMATOGRAM

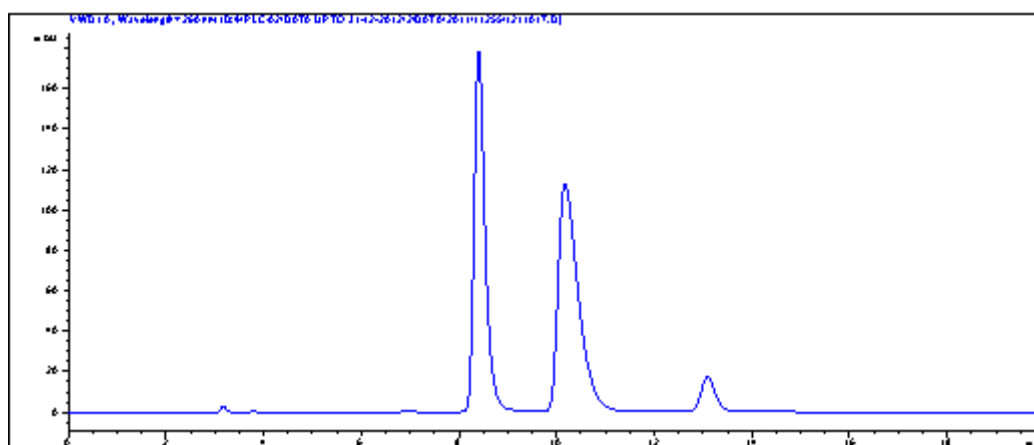

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / 1,4-dioxane / trifluoroacetic acid / diethylamine = 80 / 10 / 10 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.4	1.80	1 - 2	3.14
2	10.2	2.40		