


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Cinacalcet

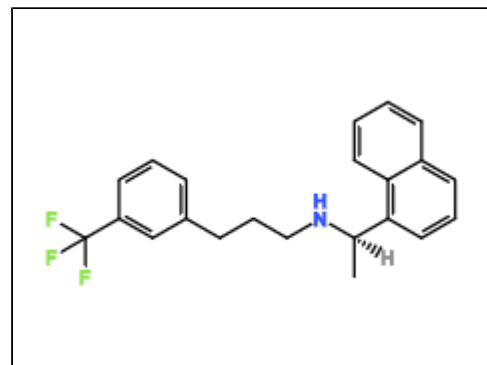

CHROMATOGRAM

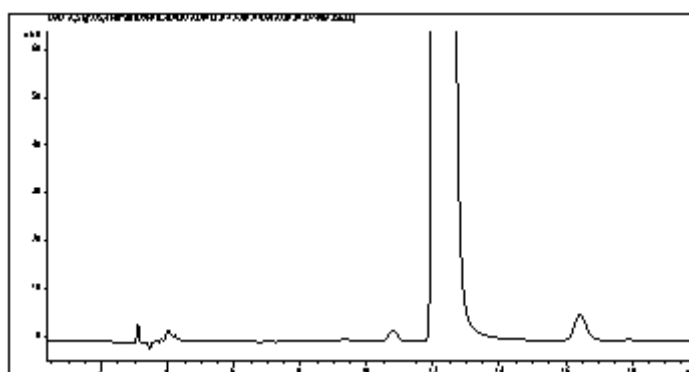

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IA-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 96 / 4 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	10.8	2.59	1 - 2	2.70
2	12.2	3.06		