


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Citalopram

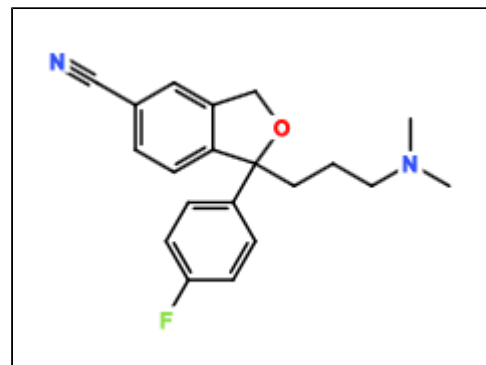

## CHROMATOGRAM

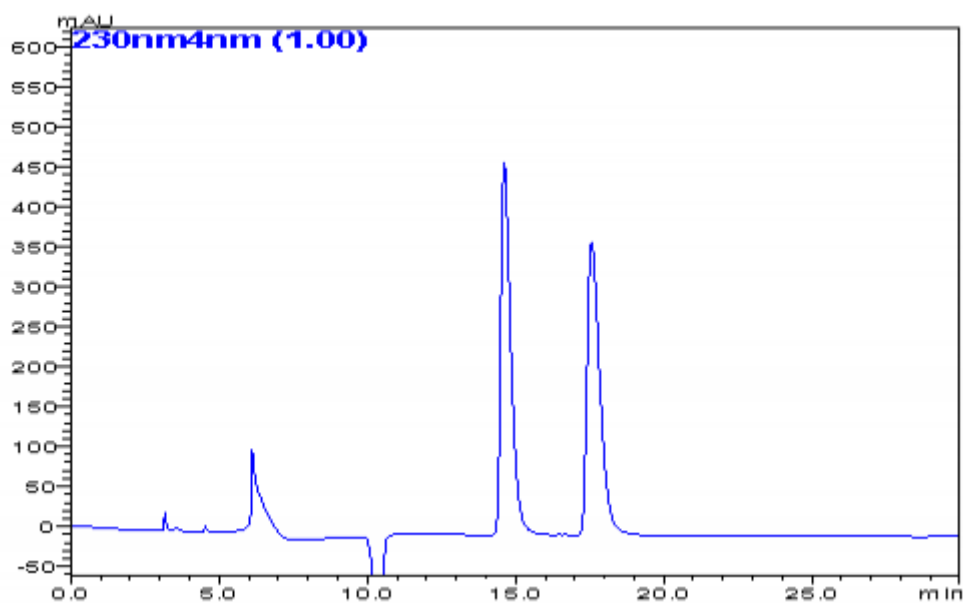

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / isopropylamine = 55 / 45 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order	API, Enantiomer		

No.	Rt	K'	Alpha Value	Resolution
1	14.6	3.87	1 - 2 1.26	3.96
2	17.6	4.85		