


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Epichlorohydrin

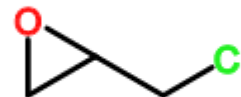

CHROMATOGRAM

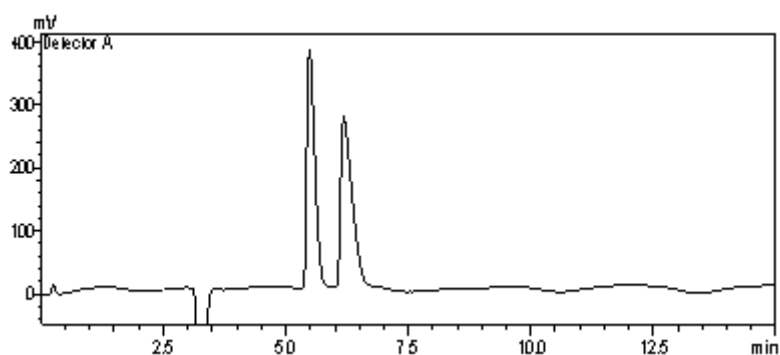

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ID 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol = 70 / 30		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.025 (mg)	Sample Concentration	
Detection	RI		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5	0.83	1 - 2	1.84
2	6.2	1.06		