


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Etomidate

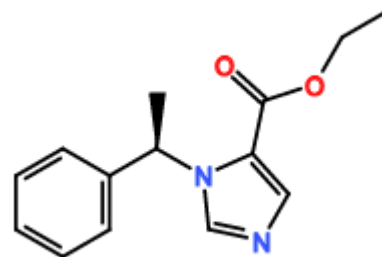

CHROMATOGRAM

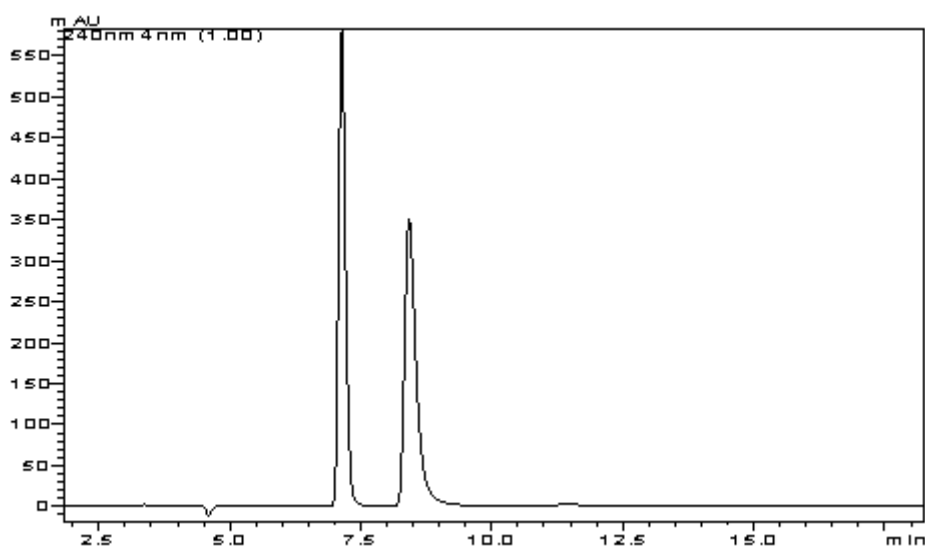

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ID 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 90 / 5 / 5 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.1	1.38	1 - 2	4.00
2	8.4	1.81		