


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Fexofenadine

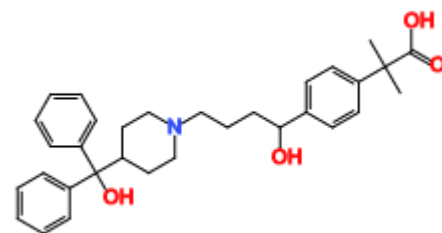

CHROMATOGRAM

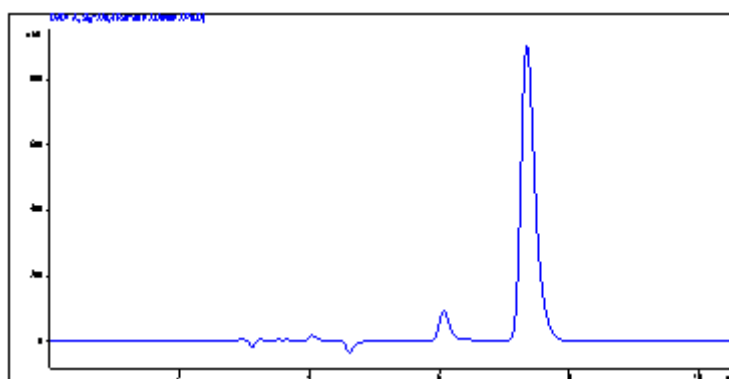

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid / diethylamine = 65 / 35 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order	Meta Isomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	6.1	1.02	1 - 2	3.79
2	7.4	1.45		