


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Flupentixol

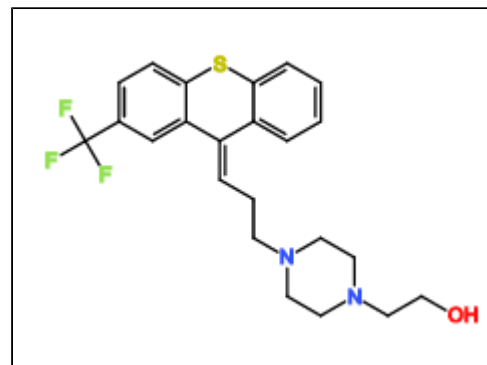

CHROMATOGRAM

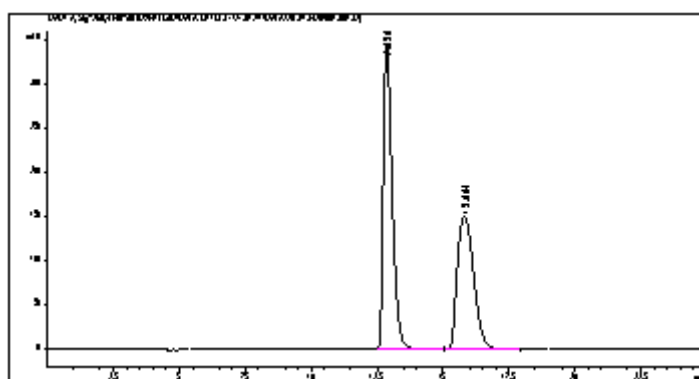

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-butanol / methanol / diethylamine = 93 / 7 / 0.5 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	15
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.9	3.29	1 - 2	1.30
2	15.8	4.27		