


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Lorcaserin

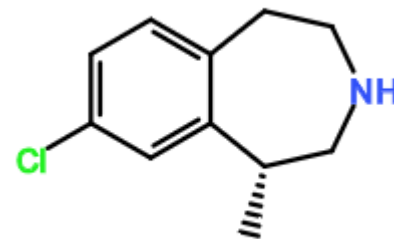

## CHROMATOGRAM

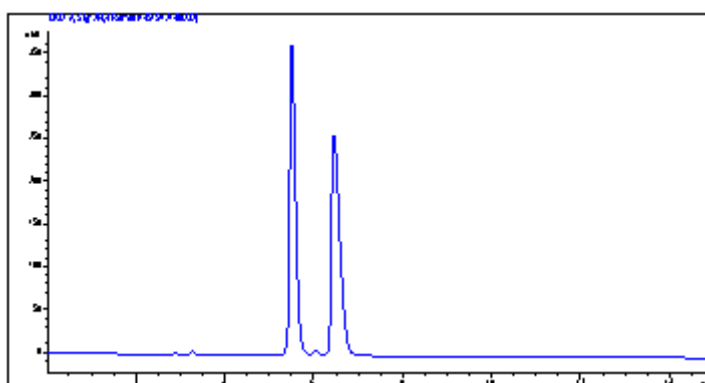

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methanol / ethylenediamine = 100 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order	S, R		

No.	Rt	K'	Alpha Value	Resolution
1	5.5	0.84	1 - 2	3.31
2	6.5	1.15		