


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Mirabegron

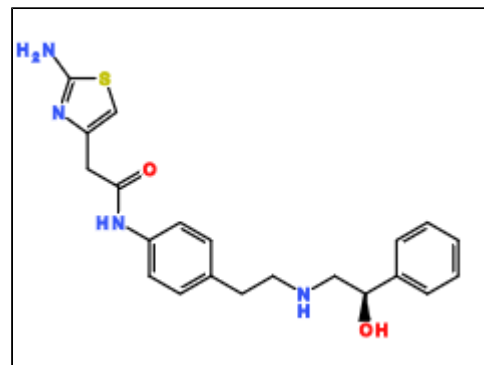

CHROMATOGRAM

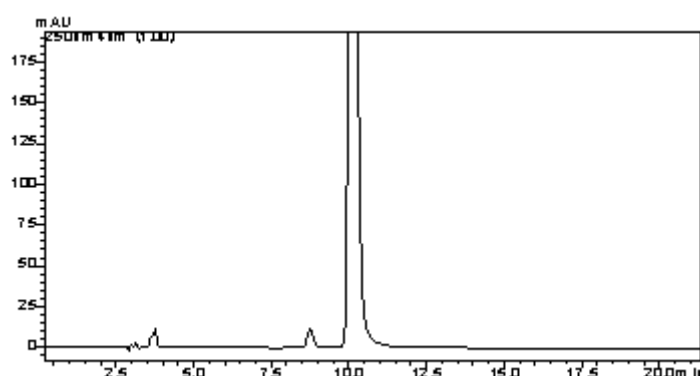

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / methanol / ethylenediamine = 30 / 30 / 40 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	8.7	1.91	1 - 2	3.98
2	10.1	2.37		