


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

O-Desmethyltramadol

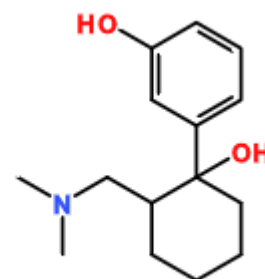

CHROMATOGRAM

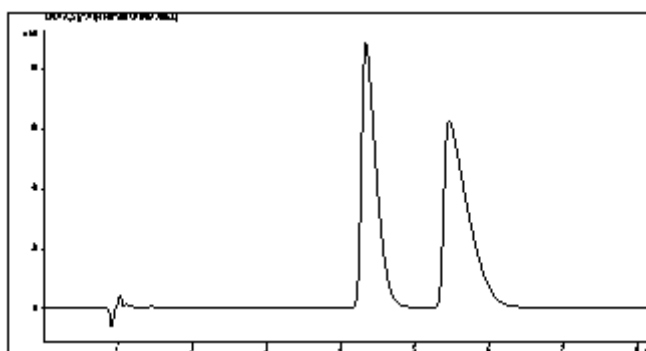

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH4OAc aq. / 2-propanol = 99 / 1		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.3	3.34	1 - 2	2.39
2	5.5	4.47		