


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Phenylalanine

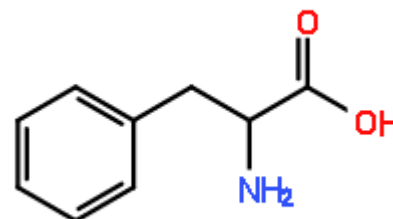

CHROMATOGRAM

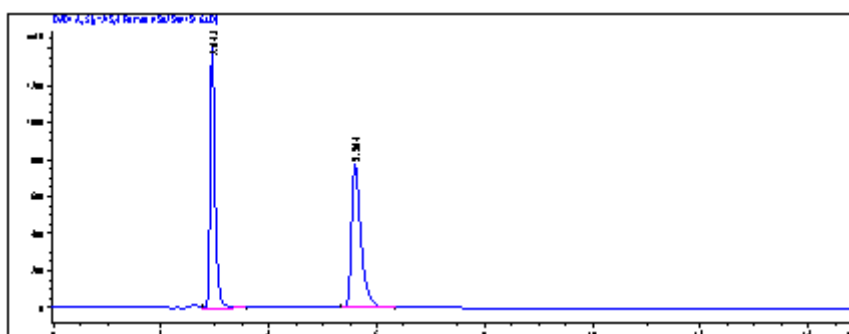

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CROWNPAK [®] CR-I(-) 3 x 150 mm (5 μm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 70 / 30		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 215 nm		
Elution Order	L, D		

No.	Rt	K'	Alpha Value	Resolution
1	2.9	0.37	1 - 2	11.80
2	5.6	1.61		