


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rebamipide

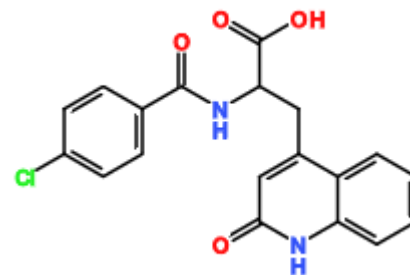

CHROMATOGRAM

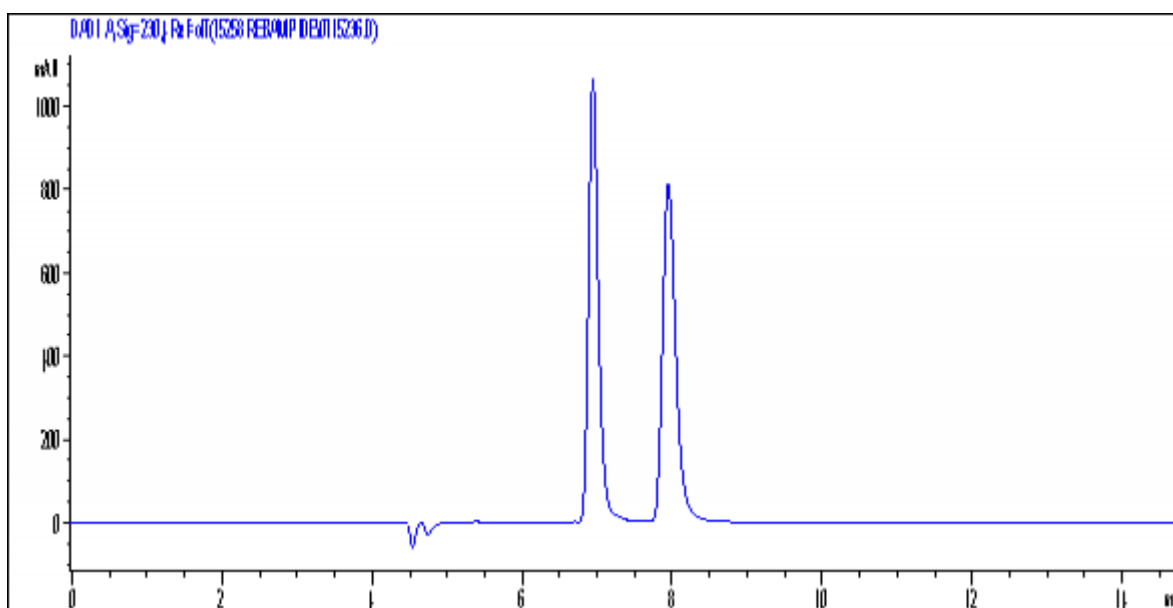

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	acetonitrile / methanol / trifluoroacetic acid = 85 / 15 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	40
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.9	0.62	1 - 2	1.38
2	8.0	0.86		