


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Saxagliptin

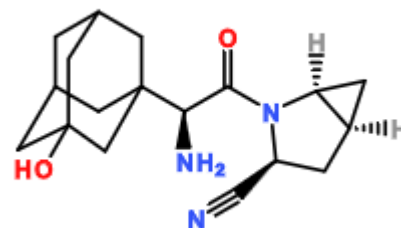

CHROMATOGRAM

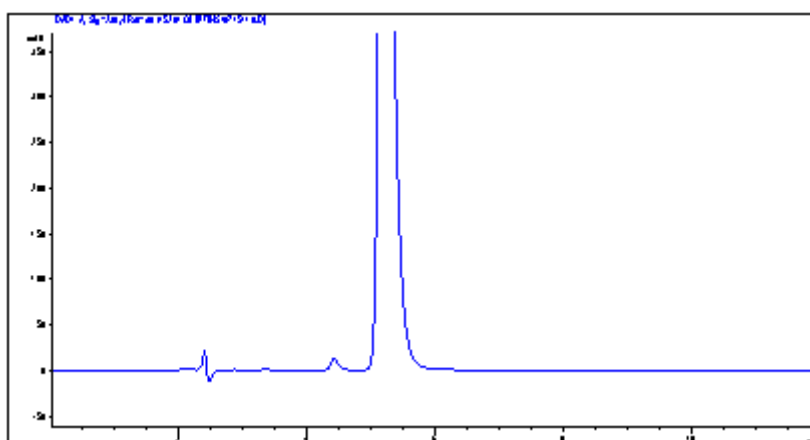

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 90 / 10		
Chromatographic Mode	Reverse Phase		

Flow rate (mL/min)	0.4	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	4.4	1.06	1 - 2	2.94
2	5.2	1.42		