


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Tryptophan

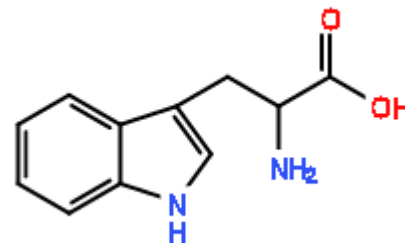

## CHROMATOGRAM

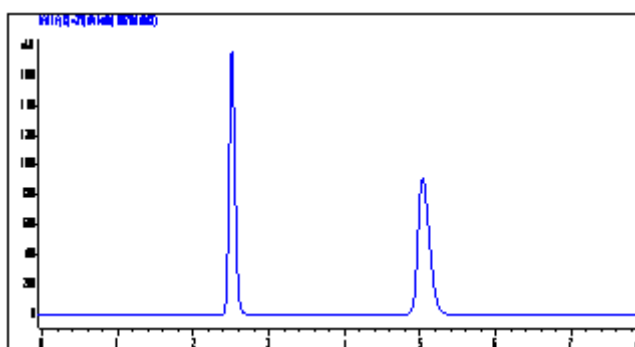

## ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.5	0.46	1 - 2	4.19
2	5.0	1.92		12.26