


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Valganciclovir

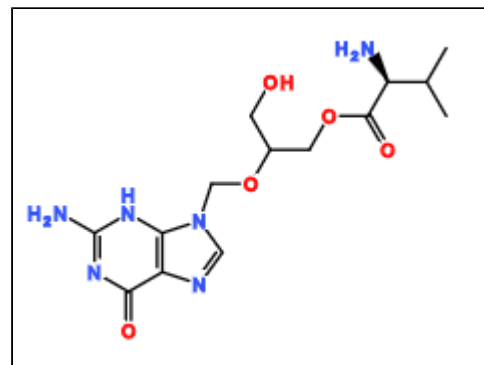

## CHROMATOGRAM

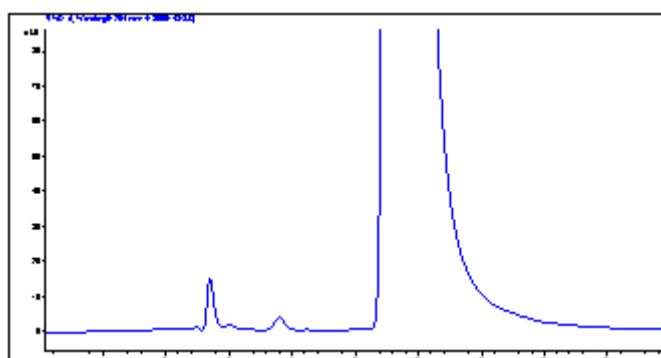

## ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 95 / 5		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order	D, L		

No.	Rt	K'	Alpha Value	Resolution
1	3.8	0.77	1 - 2	2.46
2	6.2	1.88		4.00