


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Besifloxacin

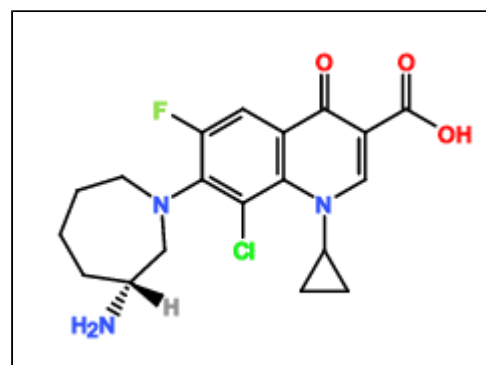

CHROMATOGRAM

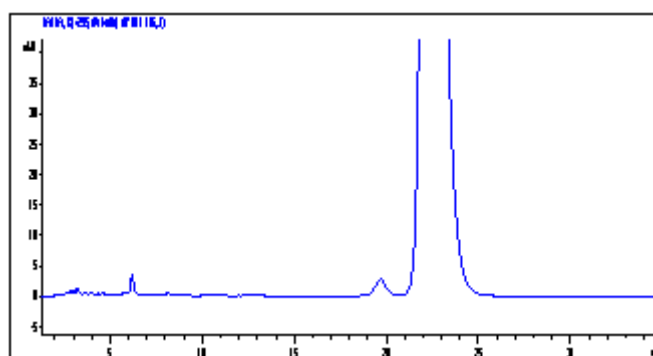

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(-)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	15
Injection Volume (µL)	0.004 (mg)	Sample Concentration	
Detection	UV-VIS 295 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	19.7	10.44	1 - 2	1.15
2	22.4	11.99		2.17