


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Dolutegravir

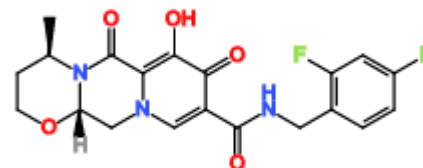

CHROMATOGRAM

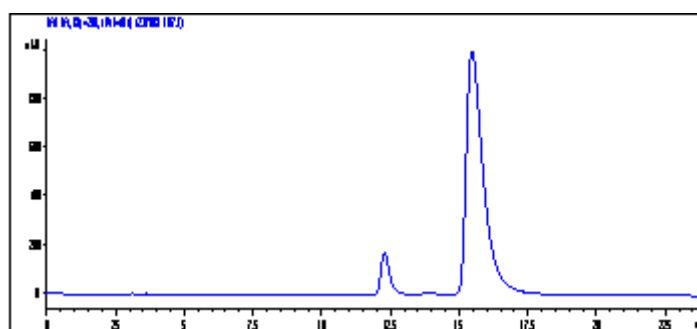

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IF-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	10mM KPF6 aq. pH2.0 / methyl tert-butyl ether / acetonitrile = 40 / 5 / 55		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	12.3	3.09	1 - 2	3.96
2	15.5	4.15		