


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Metaraminol

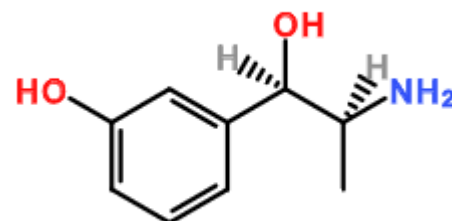

## CHROMATOGRAM

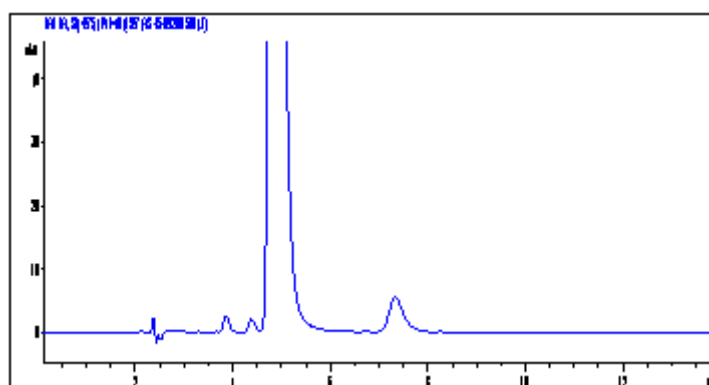

## ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	10
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order	Diastereomer, Diastereomer, API, Enantiomer		

No.	Rt	K'	Alpha Value	Resolution
1	3.9	0.80	1 - 2	1.57
2	4.4	1.03	2 - 3	2.25
3	4.8	1.25	3 - 4	5.63
4	7.3	2.41		