


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

4-Benzoyloxy-2-azetidinone

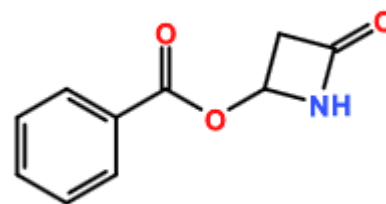

CHROMATOGRAM

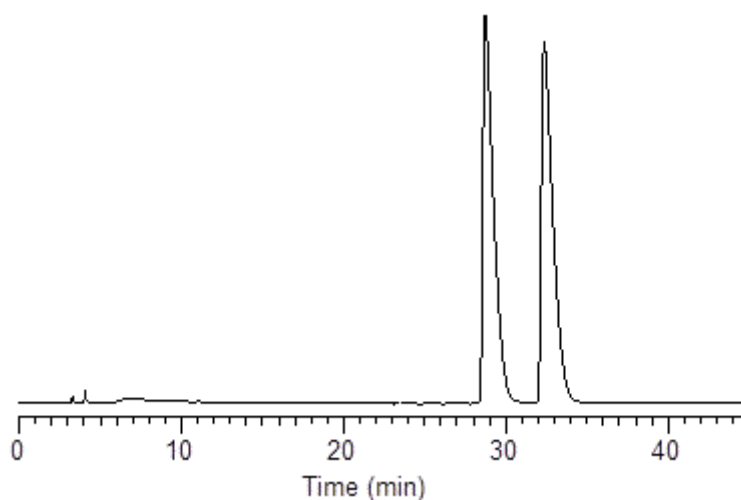

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	28.7	8.60	1 - 2	2.97
2	32.4	9.80		