


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-CBZ-Asparatic acid

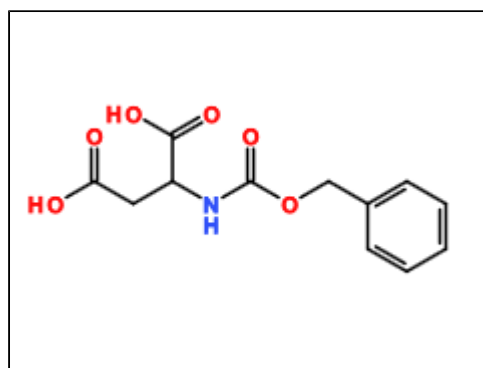

CHROMATOGRAM

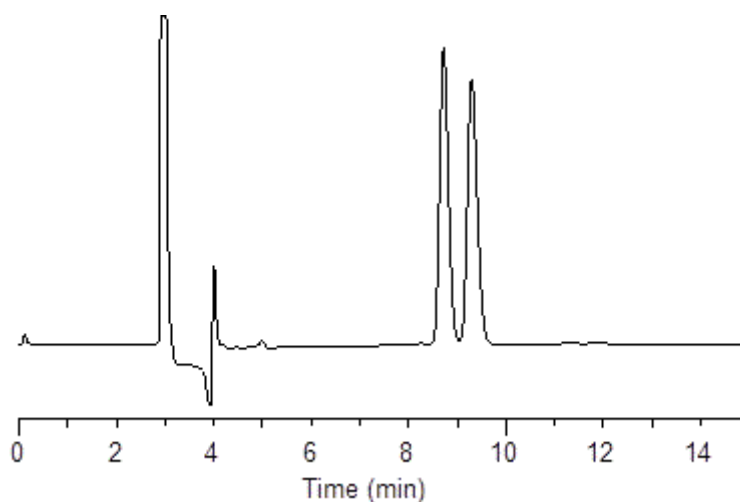

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IE 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.7	1.90	1 - 2 1.11	1.62
2	9.3	2.10		