


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Efavirenz

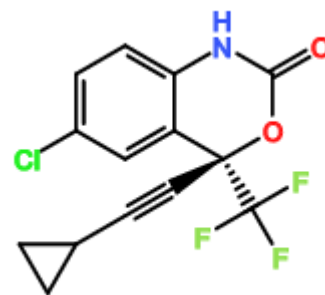

CHROMATOGRAM

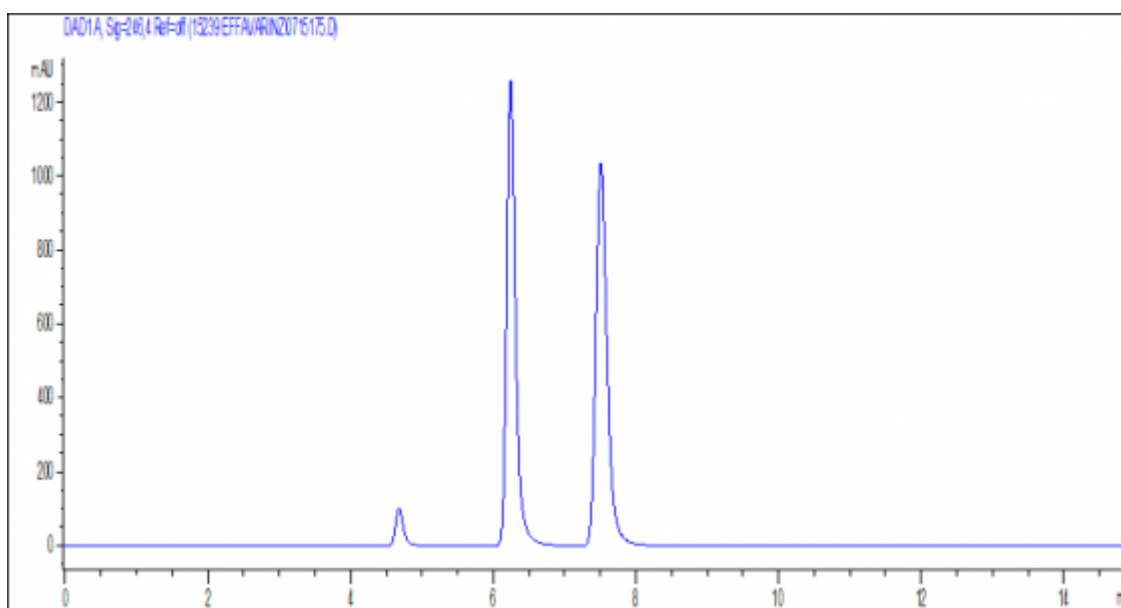

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	25mM ammonium bicarbonate aq. / methanol = 20 / 80		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 246 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.2	0.73	1 - 2	5.19
2	7.5	1.09		