


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Fosamprenavir

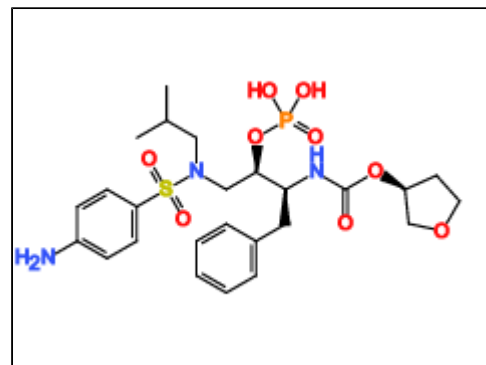

CHROMATOGRAM

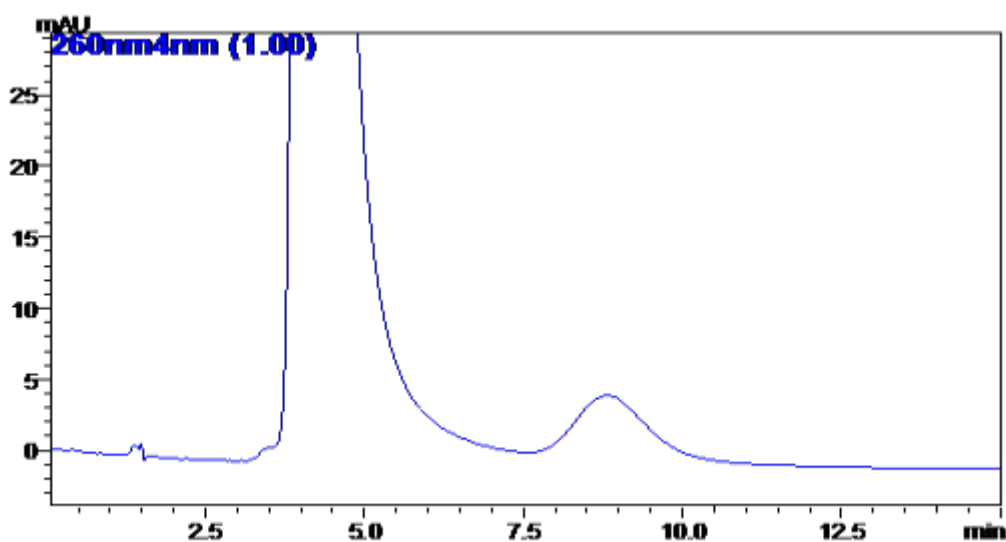

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	20mM NaH ₂ PO ₄ aq. pH6.5 / 2-propanol = 98 / 2		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.1	1.94	1 - 2	2.73
2	8.8	5.29		3.30