


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Linezolid

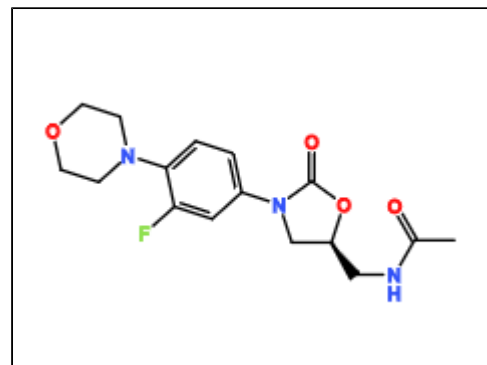

## CHROMATOGRAM

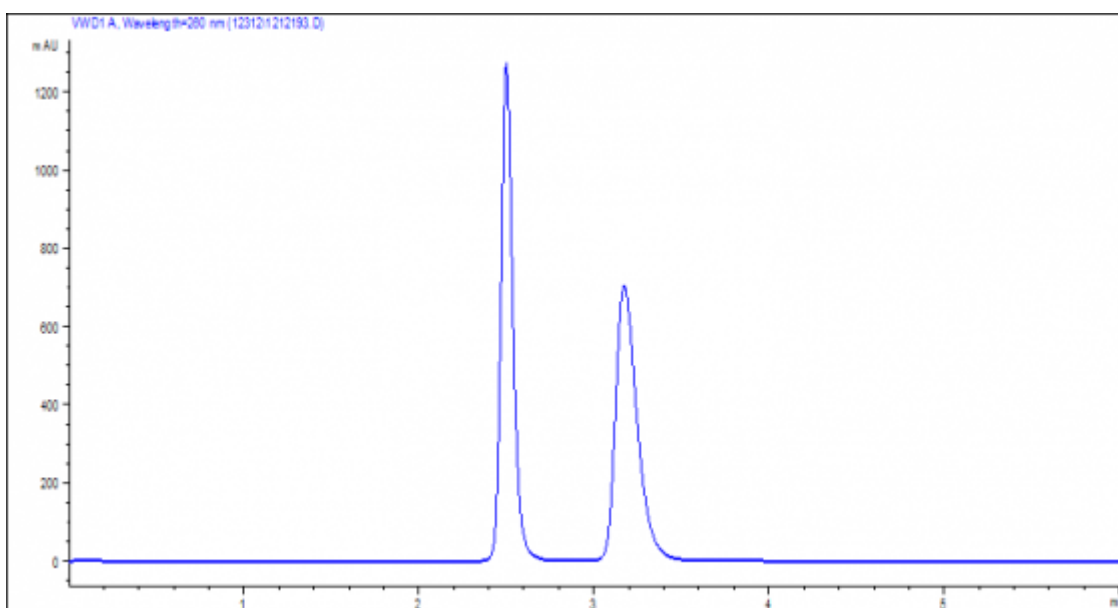

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	100mM ammonium bicarbonate aq. / acetonitrile = 5 / 95		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.5	0.37	1 - 2	3.95
2	3.1	0.71		